


Portable Electronic Irrigator Model NO.RLI502

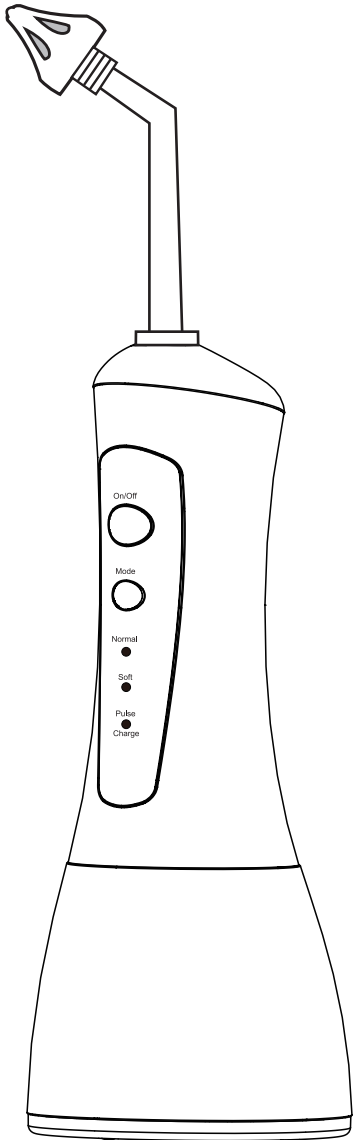

Before operating this product, please read these instructions completely and save this manual for future use.

IMPORTANT SAFEGUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

READ ALL INSTRUCTIONS BEFORE USING

DANGERS


- Toddlers, infants, people who are unable to operate should not allowed use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into tub or sink.
- Do not reach for the product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.

WARNINGS

- Electronic irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Children shall not play with the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- Regarding how to remove the battery, please see "Disposal" in page 4.
- Do not fill the tank with water warmer over 40°C. Never use the adaptor if it has a damaged cord or plug.
- Keep the cord away from heated surfaces.
- Please do not use it if you lose the wand.
- Do not submerge the device in water.
- Keep wand away from small children. The wand is small enough to be swallowed and can be become lodged in the throat.
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- See 'OtoClear Tip Instructions' for additional Warnings and Contraindications related to irrigation.
- Do not use irrigator without OtoClear Tip assembled.

SAVE THESE INSTRUCTIONS

Configuration

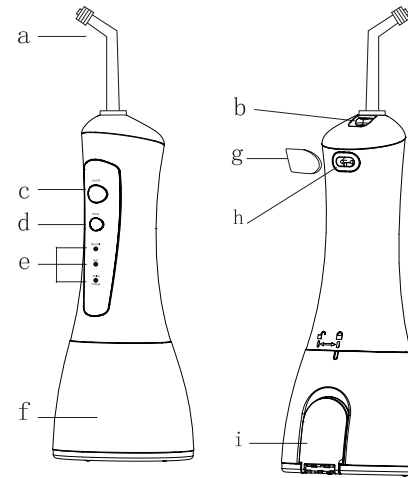

- a. Wand
- b. Wand Release Button
- c. On/Off Button
- d. Mode Button Cap
- e. Mode Indicator
- f. Water Tank
- g. Charging Base Cover
- h. Charging Base
- i. Water Tank

Note:

There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.


HOW TO USE

Charging the device (Fig. 1)

- Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).
- Insert the charging adaptor securely into the power outlet.
- The constant red LED ("Charge") indicates charging is in process.
- The constant blue LED ("Charge") indicates that the battery is fully charged.


Note:

During charging in progress, the "Charge" indicator is lit in red. After fully charged, the "Charge" indicator is lit in blue. If the "Charge" indicator is blinking red during use, it means battery exhausted, please recharge the appliance for continuous use. It takes 4 hours until fully charged. Make sure there is no liquid remaining in the electrical outlet before plugging.

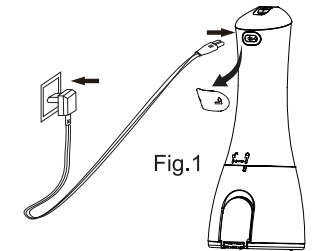


Use the water flosser


1. Open the water tank cap to fill the tank with water which temperature not exceed 40°C, then close the cap. (Fig. 2)

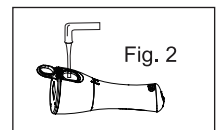


2. When first use please add small amount of water to water outlet.(Fig. 3)

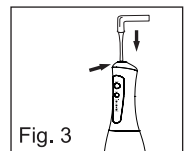

3. Attach wand onto the unit. (Fig. 4)

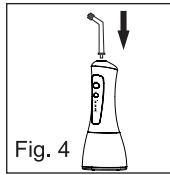


4. Press "Mode" key to select the mode you desired. During irrigating, you can change the working mode. (Fig. 5)

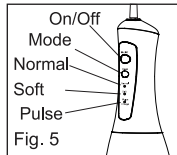


5. Hold the unit in an upright position with wand tip pointed toward yourself.
6. See the 'OtoClear Tip Instructions' for proper assembly of the OtoClear Tip onto the wand.
7. Press "On/Off" key to start irrigating.
8. Once irrigating is completed, press "On/Off" key again to turn off the device.
9. Remove and dispose of the single-use OtoClear Tip.

After using the irrigator (Fig. 6)

- Open the water tank cap and dispose of the water in the tank.
- Press "On/Off" key to turn on the unit, and drain the water off completely.
- Press "On/Off" key again to turn off the unit, and close the water tank cap.
- Wipe the device with dry cloth.

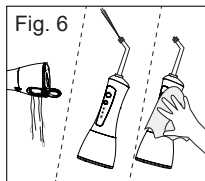

- Press the wand release button to remove the wand when you are finished. (Fig. 7)

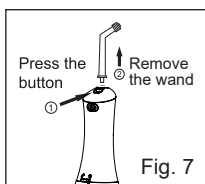

Note:

To prevent bacteria propagating, make sure there is no remaining water inside the unit. If you're not going to use it for a long time, be sure wipe and dry the appliance ASAP before store it.

Maintenance

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g: vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40°C. Keep the device away from high temperature area and avoid direct sunlight.

Main unit

- Remove and the dirt clean with a liquid soap and wipe off using a dry cloth.
- Do not submerge the main unit in water to clean.
- Dry after washing with water.

Water tank cap

Wipe off any grime and dirt with a cloth.

Wand

- Wash with water and wipe with soft cloth.
- Replace the wand every 6 months.

Water suction tube

- Wash with water and wipe with soft cloth.
- Do not bend, pull or twist the water suction tube.

Tank

- Wash with water and wipe with soft cloth.
- Dry the tank when not using for a long period.

For additional Cleaning instructions, refer to OtoClear Cleaning Instructions at bionix.com

Troubleshooting

Problem

1. It does not work
2. It only works for a few minutes even after charging
3. Water pressure is very low
4. Water does not come out

Possible Causes

- 1-1. You've just purchased the product or it has been left unused for more than 3months.
- 2-1. The battery has reached the end of its lift cycle.
- 2-2. Charging time was insufficient.
- 3-1. There is dirt stuck in the wand.
- 3-2. The wand is deformed
- 3-3. There is dirt stuck in the filter.
- 4-1. The tank is empty
- 4-2. The irrigator is tilted too much while in use
- 4-3. Valve plate slight deformation

Solution

- 1-1. Charge the unit for at least 4 hours.
- 2-1. Contact service center for support.
- 2-2. Charge the unit for at least 4 hours.
- 3-1. Replace the wand
- 3-2. Replace the wand
- 3-3. Clean the filter
- 4-1. Fill up tank
- 4-2. Use with the appliance in a vertical state
- 4-3. Add small amount of water to water outlet

